TECHNICAL DATA SHEET

THUNDERTM Total GSK3β

TR-FRET Cell Signaling Assay Kit

CATALOG NUMBERS KIT-GSK3T-100 (100 tests)

KIT-GSK3T-100 (100 tests)
KIT-GSK3T-500 (500 tests)
KIT-GSK3T-2500 (2500 tests)
KIT-GSK3T-5000 (5000 tests)
KIT-GSK3T-10000 (10000 tests)

Store at -80°C For research use only. Not for use in diagnostic procedures.



PRODUCT DESCRIPTION

This assay kit measures intracellular levels of **total GSK3** β protein in cell lysates using a simple, rapid and sensitive immunoassay based on the homogeneous (no-wash) THUNDERTM TR-FRET technology. The kit is compatible with both adherent and suspension cells.

SPECIFICITY

This assay kit contains two specific and selective antibodies that recognize **total** (both phosphorylated and unphosphorylated) $GSK3\beta$.

SPECIES REACTIVITY

Human (Swiss-Prot Acc.: P49841; Entrez-Gene Id: 2932).

Other species should be tested on a case-by-case basis.

TR-FRET ASSAY PRINCIPLE

The **Total GSK3**β assay kit is a homogeneous time-resolved Förster resonance energy transfer (TR-FRET) sandwich immunoassay (Figure 1). The THUNDER™ Cell Signaling assay workflow consists of 3 steps (Figure 2). Following cell treatment, cells are first lysed with the specific Lysis Buffer provided in the kit. Then **Total GSK3** β in the cell lysates is detected with a pair of fluorophore-labeled antibodies in a simple "add-incubate-measure" format (single-step reagent addition; no wash steps). One antibody is labeled with a donor fluorophore (Europium chelate: Eu-Abl) and the second with a farred acceptor fluorophore (FR-Ab2). The binding of the two labeled antibodies to distinct epitopes on the target protein takes place in solution and brings the two dyes into close proximity. Excitation of the donor Europium chelate molecules with a flash lamp (320 or 340 nm) or a laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which in turn emit a TR-FRET signal at 665 nm. Residual energy from the Eu chelate generates light at 615 nm. The signal at 665 nm is proportional to the concentration of **Total GSK3\beta** in the cell lysate. Data can be expressed as either the signal at 665 nm or the 665 nm/615 nm ratio.

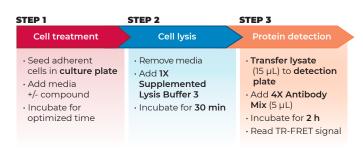


Figure 2 Assay workflow using the 2-plate (transfer) protocol.

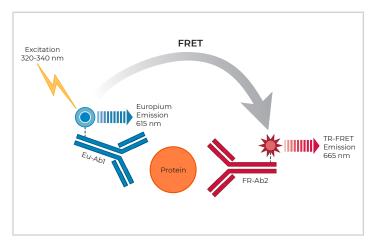


Figure 1 Schematic representation of the TR-FRET cell signaling assay principle.

KIT COMPONENTS

| | 100 points* | 500 points* |
|--|-------------|-------------|
| Eu-labeled total-GSK3β antibody (Eu-Abl) | 5 µL | 25 μL |
| Acceptor-labeled total-GSK3β antibody (FR-Ab2) | 20 μL | 100 μL |
| Lysis Buffer 3 (5X) | 1 mL | 5 mL |
| Detection Buffer (10X) | 50 μL | 250 μL |
| Positive control cell lysate | 100 μL | 200 μL |
| Phosphatase Inhibitor Cocktail (100X) | 50 μL | 250 μL |


 $^{^*}$ The number of assay points is based on an assay volume of 20 μ L in half-area 96-well or low-volume 384-well assay plates using the kit components at the recommended concentrations (refer to the User Manual).

VALIDATION DATA

This assay kit has been validated for the relative quantification of total GSK3 β in MCF7 cell lysates using the 2 plate assay protocol.


- Adherent cells were cultured overnight in a 96-well tissue culture plate (EMEM +10% FBS).
- Following cell treatment, the media was removed and cells were lysed with the 1X Lysis Buffer 3 (50 µL) supplemented with the phosphatase inhibitors sodium fluoride (1 mM) and sodium orthovanadate (2 mM).
- · Following a **30-min** incubation at room temperature (RT) on an orbital shaker (400 rpm), lysates (15 μ L) were then transferred to a 384-well assay plate followed by addition of the labeled antibodies Eu-Ab1 and FR-Ab2 (5 μ L) for detection of total GSK3 β .
- The plate was incubated at RT for **2 hours** and the TR-FRET signal was recorded at 665 and 615 nm (EnVision®; lamp excitation).

STIMULATION OF PHOSPHO-GSK3 β (S9) IN MCF7 CELLS

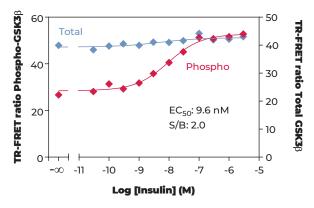

IN MCF7 CELLS


INHIBITION OF PHOSPHO-GSK3B

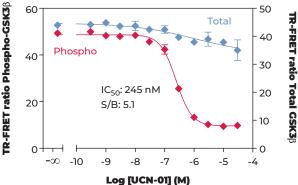

MCF7 cells (75,000 cells/well; in triplicate) were incubated with serial dilutions of Insulin for 30 min at RT. Data show that treatment of MCF7 cells with Insulin stimulates phosphorylation of GSK3 β at S9 but does not affect the levels of total GSK3 β .

MCF7 cells (75,000 cells/well; in triplicate) were incubated with serial dilutions of the UCN-01 for 2 hours at RT. Cells were then stimulated with 100 nM Insulin for 30 min at RT. Data show that treatment of MCF7 cells with UCN-01 inhibits phosphorylation of GSK3 β at S9 by Insulin, but does not have a major effect on the levels of total GSK3 β .

MCF7 CONTROL CELL LYSATE TITRATION (QC TEST)

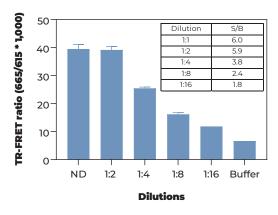

Quality Control: the Total GSK3 β assay kit is routinely tested against Insulin-treated MCF7 lysates. MCF7 cells were cultured in a T175 flask to 65% confluence and stimulated with 3 μ M Insulin for 30 min at RT. Following cell lysis using 4 mL of 1X Lysis Buffer 3, lysates were serially diluted with 1X Lysis Buffer 3 and tested in triplicate. Data show a linear relationship between lysate dilutions and TR-FRET ratio values. Note that due to the very high sensitivity of the Total GSK3 β kit, lysates from the T175 flask required at least a 1:2 pre-dilution in order to be within the dynamic assay range.



FOR MORE INFORMATION ON DEVELOPING AND OPTIMIZING TR-FRET CELL SIGNALING ASSAYS, CONSULT THE USER MANUAL.