

TECHNICAL DATA SHEET



bioauxilium

BETTER TOOLS. REAL DISCOVERIES.

THUNDER™ CHO HCP TR-FRET Biomarker Assay Kit

CATALOG NUMBERS KIT-CHOHCP-100 (100 tests)
KIT-CHOHCP-500 (500 tests)
KIT-CHOHCP-2500 (2500 tests)
KIT-CHOHCP-5000 (5000 tests)
KIT-CHOHCP-10000 (10000 tests)

Store at -80°C.
For research use only.
Not for use in diagnostic procedures.

[Visit Product Page](#)

PRODUCT DESCRIPTION

This assay kit is designed for the quantitative determination of **Host Cell Protein** from **Cultured Chinese Hamster Ovary (CHO HCP)** in cell culture supernatants or biotherapeutic protein formulations, using a simple, rapid and sensitive immunoassay based on the homogeneous (no wash) THUNDER™ TR-FRET technology. The kit does not include a **CHO HCP** standard, which must be purchased separately from Cygnus Technologies, LLC (catalogue number F018H-10).

SPECIFICITY

This assay kit contains two specific and selective antibodies that recognize **CHO HCP** in cell culture supernatant or biotherapeutic protein formulations samples.

SPECIES REACTIVITY

Chinese Hamster Ovary cells, type A (total HCP)

TR-FRET ASSAY PRINCIPLE

The THUNDER™ **CHO HCP** assay kit is a homogeneous time-resolved Förster resonance energy transfer (TR-FRET) sandwich immunoassay (Figure 1). The assay workflow consists of a one-step reagent addition (Figure 2). **CHO HCP** in the cell supernatant sample is detected with a pair of fluorophore-labeled antibodies in a simple "add-incubate-measure" format (no wash steps). One antibody is labeled with a donor fluorophore (Europium chelate; Eu-Ab1) and the second with a far-red acceptor fluorophore (FR-Ab2). The binding of the two labeled antibodies to distinct epitopes on the target cytokine takes place in solution and brings the two dyes into close proximity. Excitation of the donor Europium chelate molecules with a flash lamp (320 or 340 nm) or a laser (337 nm) triggers a FRET from the donor to the acceptor molecules, which in turn emit a TR-FRET signal at 665 nm. Residual energy from the Eu chelate generates light at 615 nm. The signal at 665 nm is proportional to the concentration of **CHO HCP** in the cell supernatant. Data can be expressed as either the signal at 665 nm or the 665 nm/615 nm ratio. A distinct advantage of the THUNDER™ TR FRET technology is that it does not require addition of potassium fluoride to ensure signal stability. The non-toxic nature of THUNDER™ makes it an ideal candidate for multiplexing with other assay technologies using a sequential assay protocol.

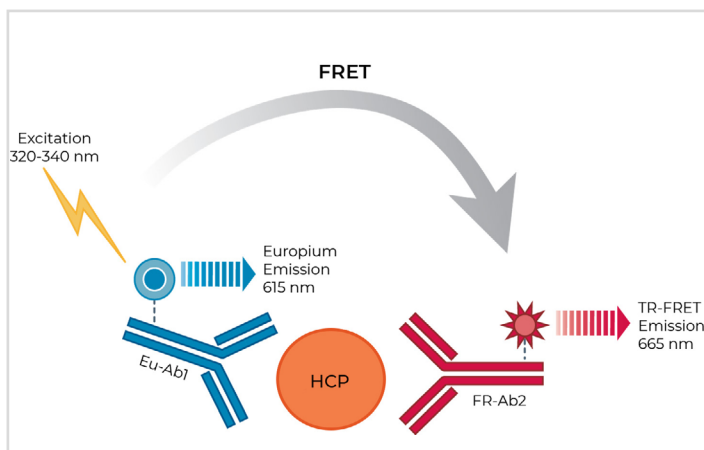


Figure 1 Schematic representation of the TR-FRET assay principle.



Figure 2 THUNDER™ CHO HCP assay workflow.

Upon receipt, store kit in the dark at -80°C.

THUNDER™ ASSAY PROTOCOL

This protocol must be read in its entirety prior to beginning the assay.

- We cannot guarantee the performance of the product outside the conditions detailed in this Assay Protocol.
- Bring all reagents to room temperature prior to running the assay.
- Centrifuge all tubes before use to improve recovery of content (2000x g, 10-15 sec).
- Use ultrapure water (Milli-Q® grade water; 18 MΩ cm) to dilute Assay Buffer.
- It is recommended to test all standards in triplicate and samples at least in duplicate.

KIT COMPONENTS

	100 points*	500 points*
Eu-labeled anti-CHO HCP antibody (Eu-Ab1)	5 µL (1 clear tube, red cap)	25 µL (1 clear tube, red cap)
Acceptor-labeled anti-CHO HCP antibody (FR-Ab2)	20 µL (1 brown tube, blue cap)	100 µL (1 brown tube, blue cap)
Assay Buffer 3 (5X)**	1 mL (1 tube, yellow cap)	3 mL (3 tubes, yellow cap)

*The number of assay points is based on an assay volume of 20 µL in low-volume 384-well assay plates using the kit components at the recommended concentrations.

**Extra Assay Buffer can be ordered separately (cat # TRFRET-AB3).

ADDITIONAL MATERIALS REQUIRED

	Recommended source	Catalog No.
CHO HCP Lysate Antigen Concentrate	Cygnus Technologies	F018H-10
Ultrapure laboratory grade water	Many options available	NA
Low-volume 384-well microplate, white	PerkinElmer Greiner Corning	6007290 784075 4513
Multi- and single-channel pipettes	Many options available	NA
Adhesive sealing film for plates	Many options available	NA
A plate reader equipped with the TR-FRET option	Many options available	NA

REAGENT PREPARATION

- The instructions described below are for testing the entire number of assay points in each kit. Adjust volumes accordingly for testing of fewer assay points.
- Prepare only as much reagent as is needed on the day of the experiment.

STEP 1 PREPARATION OF 1X ASSAY BUFFER 3

- Mix end-over-end the 5X Assay Buffer 3 before use.
- The thawed 5X Assay Buffer 3 can be stored at 4°C for 1-2 weeks. For longer periods of time, buffer should be stored at -80°C.
- Unused 1X Assay Buffer 3 can be stored at 4°C for 2 days.

For the 100-point kit

Add 1 mL of 5X Assay Buffer 3 to 4 mL of water.

For the 500-point kit

Add 3 mL of 5X Assay Buffer 3 to 12 mL of water.

STEP 2 PREPARATION OF CHO HCP STANDARD SERIAL DILUTIONS

IMPORTANT NOTE

The CHO HCP lysate is not included in the kit. The kit has been validated using CHO HCP lysate from Cygnus Technologies, LLC (Catalog number F018H-10). Other commercial CHO HCP lysates may also work but have not been tested with this kit.

- Prepare working standards just before use and use within one hour. Do not store working standards.
- Each well requires 15 µL of working standard.
- Change tip between each standard dilution.
- Mix each tube thoroughly before the next transfer.
- The following section describes the preparation of a standard curve for triplicate measurements (recommended).

NOTE

The Working Stock Solution should be used within one hour or aliquoted into screw-capped polypropylene vials and stored at -80°C for long term storage.

REAGENT PREPARATION (CONTINUED)

- A.** Create a 15 µg/mL Working Stock Solution by diluting, if necessary, the Stock Solution in either 1X Assay Buffer 3 if the standard curve is done in Assay buffer or in culture medium with 10% fetal bovine serum (FBS) if the standard curve is done in culture medium. Use this Working Standard Solution to prepare a standard dilution series.
- B.** Label 12 tubes, one for each working standard. Alternatively, serial dilutions can be conducted using a polypropylene 96-well plate.
- C.** Prepare serial dilutions for the standard curve as indicated in the table below using either 1X Assay Buffer 3 or culture medium with 10% FBS (for testing of culture supernatant samples):

Tube	Volume of CHO HCP	Vol. of 1X Assay Buffer or culture medium +10% FBS (µL)	[CHO HCP] in standard curve (g/mL in 15 µL)
1	30 µL diluted CHO HCP (1E-5 g/mL)	70	3.0E-6
2	30 µL of tube 1	60	1.0E-6
3	30 µL of tube 2	70	3.0E-7
4	30 µL of tube 3	60	1.0E-7
5	30 µL of tube 4	70	3.0E-8
6	30 µL of tube 5	60	1.0E-8
7	30 µL of tube 6	70	3.0E-9
8	30 µL of tube 7	60	1.0E-9
9	30 µL of tube 8	70	3.0E-10
10	30 µL of tube 9	60	1.0E-10
11	30 µL of tube 10	70	3.0E-11
12*	0	60	0

* In order to calculate the Limit of Detection (LOD), three additional zero standard tubes (13-15) are assayed in triplicate (12 wells in the 384-well plate).

STEP 3 PREPARATION OF 4X ANTIBODY DETECTION MIX

- Prepare and mix just before use.
- Each well requires 5 µL of 4X Antibody Detection Mix
- The unused 4X Antibody Detection Mix may be stored at 4°C for 2 days or aliquoted and stored at -80°C. Avoid repeated freeze/thaw cycles.

For the 100-point kit

NOTE

Due to the low reagent volumes in the 100-point kit, the antibodies are diluted with 1X Assay Buffer 3 directly in the vial.

- A.** Add 255 µL of 1X Assay Buffer 3 into the vial containing 5 µL of Eu Ab1 stock solution.
- B.** Add 240 µL of 1X Assay Buffer 3 into the vial containing 20 µL of FR Ab2 stock solution.
- C.** Mix gently 260 µL of pre-diluted Eu Ab1 with 260 µL of pre-diluted FR-Ab2.

For the 500-point kit

- A.** Mix gently 1,275 µL of 1X Assay Buffer 3 with the 25 µL of Eu Ab1 stock solution.
- B.** Mix gently 1,200 µL of 1X Assay Buffer 3 with the 100 µL of FR-Ab2.
- C.** Mix gently 1,300 µL of pre-diluted Eu Ab1 with 1,300 µL of pre-diluted FR Ab2.

SAMPLE PREPARATION

- Each well requires 15 µL of sample (cell supernatant).
- Just after their collection, put the samples at +4°C and assay immediately or aliquot and store -80°C. Avoid repeated freeze/thaw cycles.
- Do not use a water bath to thaw samples. Thaw samples at room temperature.
- If the analyte concentration in the sample possibly exceeds the highest point of the standard curve, prepare one or more 10-fold dilutions of the sample. Prepare serial dilutions using the same medium used to culture the cells.

ASSAY PROCEDURE

- Samples and standards must be assayed at least in duplicate each time the assay is performed.
- When loading reagents in the low-volume 384-well microplate, change tips between each standard or sample addition and after each set of reagents.
- If using a multichannel pipettor, always use a new disposable reagent reservoir.
- When reagents are added to the microplate, make sure the liquids are at the bottom of the well.
 - Add 15 µL of each working standard or 15 µL of cell supernatant sample.
 - Add 5 µL of 4X Antibody Detection Mix (Eu-Ab1 + FR-Ab2) to each of the assay wells.
 - Cover the plate with a plate sealer and incubate 18 hours at room temperature.
 - Gently remove the adhesive plate sealer. Read plate on a TR-FRET compatible microplate reader.

NOTE *The same plate can be read several times without any negative effect on the assay performance.*

TR-FRET PLATE READER SETTINGS

For optimal assay performance, we recommend that you read the TR FRET assays at two wavelengths, detecting both the emission from the Europium chelate donor fluorophore at 615 nm, and the acceptor fluorophore at 665 nm. The following instrument settings are provided as guidelines.

Parameter	TR-FRET Compatible Plate Reader*	
	Flash lamp excitation	Laser excitation
Excitation filter	320 nm (or 340 nm)	Not applicable
Emission filter	615 nm (or 620 nm)	615 nm (or 620 nm)
Delay time	90 µs	50 µs
Flash energy level	100% or High	100%
Number of flashes	100-200	20
Window (integration time)	300 µs	100 µs

* These settings are provided as guidelines only. Settings should be optimized for each reader. For optimal settings, visit bioauxilium.com/resources.

DATA ANALYSIS

- TR-FRET data are typically calculated and presented ratiometrically using the following formula: $[(665 \text{ nm}/615 \text{ nm}) \times 1,000]$
- Alternatively, the signals at 665 nm can be used directly to process your data.
- Calculate the TR-FRET ratio for each well.
- Since TR-FRET signal is read in a time-resolved mode, background subtraction is usually not necessary.
- Create a standard curve by analyzing data according to a nonlinear regression using the 4 parameter logistic equation (sigmoidal dose-response curve with variable slope) and a $1/Y^2$ data weighting. For correct analysis, the highest standard point should be removed from the curve if it is present after the hook point (i.e., if it shows lower counts compared to the next standard).
- Do not use a value of zero for the zero standard. Use instead the actual background values obtained with the zero standard.
- If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

TYPICAL STANDARD CURVE

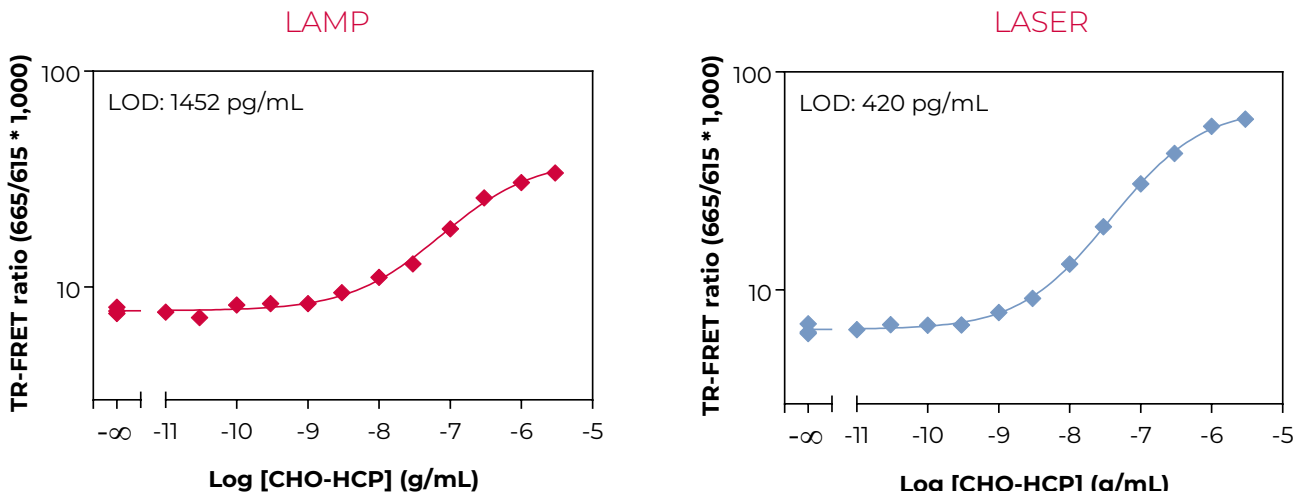


Figure 3 Typical nonlinear CHO HCP standard curves in Assay Buffer 3. Left: lamp excitation, Right: laser excitation. Data represent the mean ± standard deviation of triplicate measurements (3 wells) for each standard. The data was generated using a low-volume 384-well white plate read on the EnVision® equipped with TR-FRET option.

VALIDATION DATA

The THUNDER™ CHO HCP kit has been validated using CHO HCP lysate antigen concentrate from Cygnus Technologies, LLC (Catalog number F018H-10; NOT included in the kit).

CHO HCP THUNDER™ ASSAY SUMMARY		
Assay Type	Homogeneous sandwich immunoassay	
Format	384-well plate	
Hands-on time	15 minutes	
Incubation time	18 hours	
Sample type (volume)	Cell culture supernatant (15 µL)	
Specificity	CHO HCP, type A (total HCP)	
	Laser excitation	Lamp excitation
Limit of Detection (LOD)	420 pg/mL	1 452 pg/mL
Lower Limit of Quantification (LLOQ)	3 048 pg/mL	12 524 pg/mL
EC ₅₀	144 ng/mL	141 ng/mL
Dynamic Range	420 – 30,000,000 pg/mL	1452 – 30,000,000 pg/mL

SENSITIVITY

The Limit of Detection (LOD) was calculated by adding two standard deviations to the mean background counts of 12 zero standard replicates and calculating the corresponding concentration on the standard curve.

The Lowest Level of Quantitation (LLOQ) was calculated by adding 10 standard deviations to the mean background counts of 12 zero standard replicates and calculating the corresponding concentration on the standard curve.

	ASSAY BUFFER 3		F12K +10% FBS		RPMI +10% FBS	
	Laser	Lamp	Laser	Lamp	Laser	Lamp
# of standard curves	10	10	3	3	3	3
LOD (pg/mL)	420	1 452	2 353	2 252	837	1 302
LLOQ (pg/mL)	3 048	12 524	15 426	33 216	4 779	13 033

PRECISION

INTRA-ASSAY PRECISION

Three samples of known concentration (low, medium, high) were tested eighteen times each on the same plate to assess intra-assay precision. The assay was run in a 384-well plate with standards and samples prepared in Assay Buffer 3. The resulting low %CV (<10%) for each sample indicates good reproducibility within an assay.

Sample concentration	Low	Med	High
# of samples	18	18	18
Mean (pg/mL)	23 162	77 008	250 573
SD (pg/mL)	1 857	5 175	21 839
CV (%) Intra-Assay	8.0	6.7	8.7

INTER-ASSAY PRECISION

Three samples of known concentration were tested in separate assays to assess inter-assay precision. A total of eight independent experiments were performed by two operators using two different kit lots: The assay was run in 384-well plates with standards and samples prepared in Assay Buffer 3 and tested in triplicate. The resulting low %CV (<15%) for each sample indicates good reproducibility between assays.

Sample concentration	Low	Med	High
# of samples	24	24	24
Mean (pg/mL)	25 831	86 513	287 334
SD (pg/mL)	3 798	8 882	30 460
CV (%) Intra-Assay	14.7	10.3	10.6

RECOVERY

To assess the accuracy of the assay, spike-to-recovery experiments were conducted. Each assay consisted of one standard curve and three spike concentrations (low, medium and high), assayed in triplicate in Assay Buffer 3, F12K and RPMI (supplemented with 10% FBS). The average recovery was calculated from the average of 3 wells and compared to a standard curve prepared in the corresponding diluent. Data show recovery values within the range of 80–120% demonstrating assay accuracy.

	Spike level	Expected (pg/mL)	Observed (pg/mL)	Average % Recovery
Assay Buffer 3	High	300 000	358 995	120
	Medium	100 000	111 276	111
	Low	30 000	28 227	94
F12K + 10% FBS	High	300 000	286 226	95
	Medium	100 000	80 040	80
	Low	30 000	27 121	90
RPMI + 10% FBS	High	300 000	321 487	107
	Medium	100 000	95 146	95
	Low	30 000	28 486	95

LINEARITY

To assess the linearity of the assay, samples were spiked with a high concentration of analyte in culture medium supplemented with 10% FBS and diluted with the same culture medium with 10% FBS. The assay was run in a 384-well plate with standards prepared in culture media with 10% FBS. Each sample was measured in triplicate. Data show recovery values within 80-122% confirming assay linearity.

	Expected (pg/mL)	Observed (pg/mL)	Average % Recovery
F12K + 10% FBS	1 500 000	1 549 706	103
	750 000	856 148	114
	375 000	429 618	115
	187 500	206 402	110
	93 750	112 777	120
	46 875	54 129	115
	23 438	23 545	100
	11 719	11 361	97
	5 859	5 048	86
RPMI + 10% FBS	1 500 000	1 321 038	88
	750 000	782 450	104
	375 000	328 589	88
	187 500	174 417	93
	93 750	84 104	90
	46 875	42 700	91
	23 438	21 593	92
	11 719	9 447	81
	5 859	4 774	81
	2 930	2 943	100
	1 465	1 740	119

